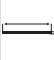
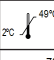
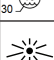
Accoat


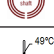
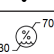
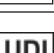
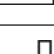
Guide Wire

EN	<i>English</i>	Guide Wire
BG	<i>Bulgarian/Български език</i>	Водач
ZH	<i>Chinese/简体中文</i>	导丝
HR	<i>Croatian/Hrvatski</i>	Žica vodilica
CS	<i>Czech/Česky</i>	Vodící drát
DA	<i>Danish/Dansk</i>	Guide Wire
NL	<i>Dutch/Nederlands</i>	Voerdraad
ET	<i>Estonian/Eesti</i>	Juhtetraat
FI	<i>Finnish/Suomi</i>	Ohjainvaijeri
FR	<i>French/Français</i>	Fil-guide
DE	<i>German/Deutsch</i>	Führungsdraht
EL	<i>Greek/Ελληνικά</i>	Συρμάτινος οδηγός
HU	<i>Hungarian/Magyar</i>	Vezetődrót
IT	<i>Italian/Italiano</i>	Filo guida
LV	<i>Latvian/Latviešu</i>	Vadītājstīga
LT	<i>Lithuanian/Lietuvių k.</i>	Kreipiamoji viela
NO	<i>Norwegian/Norsk</i>	Ledesonde
PL	<i>Polish/Polski</i>	Prowadnik
PT	<i>Portuguese/Potuguês</i>	Fio-guia
RO	<i>Romanian/Română</i>	Fir de ghidaj
RU	<i>Russian/Русский</i>	Проволочный проводник катетера
SL	<i>Slovenian/Slovenščina</i>	Vodilna žica
ES	<i>Spanish/Español</i>	Alambre guía
SV	<i>Swedish/Svensk</i>	ledare
TR	<i>Turkish/Türkçe</i>	Kılavuz Tel

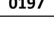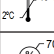
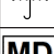
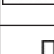
EN: Explanation of symbols used on the package labels - BG: Значение на символите върху етикетите на опаковката - ZH: 包装标签上的标志说明 - HR: Objašnjenja simbola koji se koriste na naljepnici na pakiranju - CS: Vysvětlení symbolů použitých na štítcích obalů - DA: Forklaring af symboler anvendt på emballagens etiketter

	EN: Caution - BG: Внимание - ZH: 注意 - HR: Oprez - CS: Upozornění - DA: Forsigtig
	EN: Consult instructions for use or consult electronic instruction for use - BG: Направете справка с инструкциите за употреба или с електронните инструкции за употреба - ZH: 参阅使用说明或参阅电子版使用说明 - HR: Pogledajte upute za uporabu u tiskanom ili elektroničkom obliku - CS: Přečtěte si návod k použití nebo elektronický návod k použití - DA: Se brugsanvisningen, eller se den elektroniske brugsanvisning
	EN: Contents - BG: Съдържание - ZH: 内含物 - HR: Sadržaj - CS: Obsah - DA: Indhold
	EN: Guide wire length and diameter - BG: Дължина и диаметър на водача - ZH: 导丝长度及直径 - HR: Duljina i promjer žice vodilice - CS: Délka a průměr vodičoho drátu - DA: Guidewirens længde og diameter
	EN: Tip shape: Straight - BG: Форма на върха: Прав - ZH: 末端形状：直形 - HR: Oblik vrha: Ravan - CS: Tvar hrotu: Rovný - DA: Spidsform: Lige
	EN: Tip shape: J - BG: Форма на върха: J-образен - ZH: 末端形状：J形 - HR: Oblik vrha: J - CS: Tvar hrotu: Zakřivený - DA: Spidsform: J-format
	EN: Tip shape: J and straight - BG: Форма на върха: J-образен и прав - ZH: 末端形状：J形和直形 - HR: Oblik vrha: J i ravan - CS: Tvar hrotu: Zakřivený a rovný - DA: Spidsform: J-format og lige
	EN: Shaft stiffness: Standard - BG: Твърдост на стеблото: стандартно - ZH: 轴刚度：标准 - HR: Čvrstoća osovine: Standardna - CS: Tuhost dříku: Standardní - DA: Skafkets stivhed: Standard
	EN: Shaft stiffness: Stiff - BG: Твърдост на стеблото: твърдо - ZH: 轴刚度：硬性 - HR: Čvrstoća osovine: Čvrsta - CS: Tuhost dříku: Tuhý - DA: Skafkets stivhed: Stiv
	EN: Temperature limit - BG: Диапазон за температурата - ZH: 温度限制 - HR: Ograničenje temperature - CS: Teplotní mez - DA: Temperaturgrænse
	EN: Humidity limit - BG: Диапазон за влажността - ZH: 湿度限制 - HR: Ograničenje vlažnosti - CS: Mez vlhkosti - DA: Fugtighedsgrænse
	EN: Keep away from sunlight - BG: Да се пази от слънчева светлина - ZH: 避免日光照射 - HR: Držati podalje od sunčeve svjetlosti - CS: Chraňte před slunečním světlem - DA: Holdes væk fra sollys
	EN: Keep dry - BG: Да се съхранява на сухо място - ZH: 保持干燥 - HR: Proizvod se mora održavati suhim - CS: Uchovávejte v suchu - DA: Opbevares tørt
	EN: Medical Device - BG: Медицинско изделие - ZH: 医疗器械 - HR: Medicinski uređaj - CS: Zdravotnický prostředek - DA: Medicinsk udstyr
	EN: Sterilized by ethylene oxide gas. Single sterile barrier system with protective packaging inside. - BG: Стерилизирано с газ етиленов оксид. Единична стерилна бариерна система с предпазна опаковка вътре. - ZH: 用环氧乙烷气体消毒。内部带有保护性包装的单一无菌屏障系统。 - HR: Sterilizirano etilen oksidom. Sustav s jednom sterilnom barijerom i unutarnjim zaštitnim pakiranjem. - CS: Sterilizováno plynným etylenoxidem. Systém jedné sterilní bariéry s ochranným obalem uvnitř. - DA: Steriliseret med ætylenoxid. Sterilt enkeltbarrieresystem med beskyttende emballage indeni.
	EN: Do not use if pack is open or damaged - BG: Да не се използва, ако опаковката е отворена или повредена - ZH: 如已开封或损坏，请勿使用。 - HR: Ne koristite ako je pakiranje otvoreno ili oštećeno - CS: Nepoužívejte, pokud je obal otevřený nebo poškozený - DA: Må ikke anvendes, hvis emballagen er åbnet eller beskadiget
	EN: Use by date - BG: Срок на годност - ZH: 使用日期 - HR: Krajnji rok uporabe - CS: Datum použitelnosti - DA: Sidste anvendelsesdato
	EN: Do not re-use - BG: Да не се използва повторно - ZH: 请勿重复使用 - HR: Nije namijenjeno za ponovnu uporabu - CS: Nepoužívejte opakovaně - DA: Må ikke genbruges
	EN: Do not re-sterilize - BG: Да не се стерилизира повторно - ZH: 请勿重新消毒 - HR: Nemojte ponovno sterilizirati - CS: Znovu nesterilizujte - DA: Må ikke gensteriliseres
	EN: Non-Pyrogenic - BG: Непириогенно - ZH: 无热原 - HR: Nepirogeno - CS: Nepyrogenní - DA: Ikke-pyrogen
	EN: Reference number - BG: Референтен номер - ZH: 产品编号 - HR: Referentni broj - CS: Referenční číslo - DA: Referencenummer
	EN: Lot number - BG: Партиден номер - ZH: 产品批号 - HR: Broj serije - CS: Číslo šarže - DA: Partinummer
	EN: UDI - Unique Device Identification - BG: UDI – уникален идентификационен номер на изделието - ZH: UDI – 唯一器械标识 - HR: UDI – Jedinstveni identifikator uređaja - CS: UDI – Jedinečná identifikace prostředku - DA: UDI – unik udstyrsidentifikation
	EN: Country of manufacturer (PL=Poland) and manufacturing date - BG: Страна на производство (PL = Poland – Полша) и дата на производство - ZH: 制造商国家 (PL = 波兰) 和制造日期 - HR: Država (PL = Poljska) i datum proizvodnje - CS: Země výrobce (PL=Polsko) a datum výroby - DA: Fremstillingsland (PL=Polen) og fremstillingsdato
	EN: Manufacturer - BG: Производител - ZH: 制造商 - HR: Proizvođač - CS: Výrobce - DA: Producent
	EN: CE mark - BG: CE маркировка - ZH: CE 认证标志 - HR: Oznaka CE - CS: Značka CE - DA: CE-mærke

NL: Verklaring van de symbolen op de verpakkingsetiketten - ET: Pakendite etiketid el kasutatavate sümbolite selgitused - FI: Pakkauksen etiketeissä käytettyjen symbolien selitykset - FR: Explication des symboles utilisés sur les étiquettes de l'emballage - DE: Erläuterung der auf der Verpackung verwendeten Symbole

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Accoat Guide Wire

Description:

The device is an EO sterilized and non-pyrogenic guide wire.

Product name	Subtype	Tip subtype	Available variants
Accoat Guide Wire	Seldinger	<ul style="list-style-type: none"> • Straight • J-tip • Double ended 	<ul style="list-style-type: none"> • Stainless steel with 2 coating variants: with/ without PTFE • Guide wire diameter • Tip shape • Guide wire stiffness • Guide wire length
	Seldinger Movable	<ul style="list-style-type: none"> • Straight • J-tip 	

For further information, see SP Medical's website.

Intended purpose:

The Accoat Guide Wire is indicated for general intravascular and coronary arterial use to aid in the selective placement of interventional devices during diagnostic and / or therapeutic procedures.

Contraindications:

Not for use in the cerebral vasculature.

It is always the physician's responsibility to determine and ensure the patient's suitability for the procedure where the Accoat Guide Wire is used.

Warnings:

1. SINGLE USE. This device is intended for single use only. Reusing the device involves high risk of contamination and locking of the guide wire inside the interventional device or human body due to guide wire kinking.
2. This device is provided in sterile condition. Do not re-sterilize or re-use the device.
3. Do not withdraw PTFE coated guide wire through a metal cannula.
4. Never advance or withdraw the guide wire against resistance until the cause of the resistance is determined by fluoroscopy.
5. Do not attempt to move the guide wire without observing the resultant tip response.
6. To be used before the expiry date stated on the package.

Precautions:

1. Prior to opening, the sterile package should be carefully inspected to see if it is still intact. The inspection must be carried out just before the medical device is used and in good light conditions. Sealing quality on the transparent side of the peel pouch must be carefully inspected for channels and cracks. Both front and back side are to be inspected for pinholes or other packaging defects such as abrasions or cracks. Use a magnification tool if in doubt of sterile barrier integrity.
2. Do not use if the package / sterile barrier is broken / damaged.
3. Prior to use, carefully inspect the guide wire for bends, kinks or other damages.
4. Do not use damaged product.
5. The device should only be used by experienced physicians, trained in invasive techniques, the use of guide wires and familiar with side effects and hazards commonly associated with interventional procedures.
6. When using movable-core guide wires, do not advance the movable core when the tip is in a curved shape.
7. Never twist or apply excessive force as the core might penetrate the coil and cause vessel damage.
8. To avoid arterial damage, make sure that the movable core is sufficiently retracted when the device is in the proximity of the heart.
9. The Accoat Guide Wire contains a metallic core, do not use with any inappropriate equipment (e.g., MRI).

Potential side effects:

Possible complications include, but are not limited to the following:

- Vessel / tract trauma
- Vessel / tract narrowing / occlusion
- Infection and inflammation
- Allergic reaction

Serious incident:

Every serious incident that has occurred in relation to the device should be reported without delay to the manufacturer and the competent authority of the Member State in which the user and / or patient is established. For further information please contact the manufacturer.

Compatibility:

The Accoat Guide Wire may be used in combination with interventional devices in CathLab environment. The end-user is solely responsible for selection of appropriate device. There should be at least 0.0004" (0,01 mm) clearance between the lumen of the interventional device and the guide wire.

Preparations for use:

1. Carefully remove the guide wire from the dispenser.
2. Inspect the guide wire thoroughly to make certain that it is not kinked or otherwise damaged.
3. Confirm the compatibility of the guide wire diameter with the interventional device before actual use.

Directions for use:

1. When introducing the guide wire into the interventional device ensure that at least 2 centimeters of guide wire extend from the proximal hub. This will prevent the guide wire from slipping inside e.g., the catheter.
2. To aid in the selective placement of the interventional device into a particular vessel, gently rotate the proximal end of the guide wire as it is advanced forward.
3. To prevent contrast agent crystallization / clotting, a continuous saline flush should be maintained between the interventional device and the guide wire during the procedure. The size of the syringe used to flush the lumen should be adapted to the length and diameter of the interventional device.
4. Between uses, during the same procedure, place the guide wire in a saline-filled container, or fill the dispenser with saline, using the luer lock provided in the package and reinsert the guide wire in the dispenser, distal end first. Make sure to leave a segment of the proximal end outside the dispenser.

Storage:

To be stored under cool, dark and dry conditions. To be kept away from heat. Humidity limits: 30% - 70%. Temperature limits: 2°C - 49°C

Expiry date:

The Accoat Guide Wire must be used no later than stated on the label.

Disposal:

Disposal shall be handled in accordance with accepted medical practice and applicable local laws and regulations for hazardous waste.

Safety information:

The IFU is delivered in hard copy together with products and kept updated and available for electronic download on SP Medical's website. Summary of safety and clinical performance can be obtained via Eudamed using SP Medical as legal manufacturer and the basic UDI-DI: 5710086-GWAC. <https://ec.europa.eu/tools/eudamed/#/screen/search-device>

Disclaimer of liability:

SP Medical A/S is not liable for defects / deterioration resulting from abnormal use or modifications made to the product and under these circumstances not covered by the guarantee. SP Medical A/S disclaims liability for direct or indirect injuries that may occur as a consequence of the product being modified or wrongly used.

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