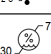
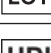
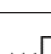
# Accoat

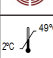
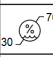
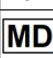
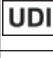
## Guide Wire

EN	<i>English</i>	<b>Guide Wire</b>
BG	<i>Bulgarian/Български език</i>	<b>Водач</b>
ZH	<i>Chinese/简体中文</i>	<b>导丝</b>
HR	<i>Croatian/Hrvatski</i>	<b>Žica vodilica</b>
CS	<i>Czech/Česky</i>	<b>Vodící drát</b>
DA	<i>Danish/Dansk</i>	<b>Guide Wire</b>
NL	<i>Dutch/Nederlands</i>	<b>Voerdraad</b>
ET	<i>Estonian/Eesti</i>	<b>Juhtetraat</b>
FI	<i>Finnish/Suomi</i>	<b>Ohjainvaijeri</b>
FR	<i>French/Français</i>	<b>Fil-guide</b>
DE	<i>German/Deutsch</i>	<b>Führungsdraht</b>
EL	<i>Greek/Ελληνικά</i>	<b>Συρμάτινος οδηγός</b>
HU	<i>Hungarian/Magyar</i>	<b>Vezetődrót</b>
IT	<i>Italian/Italiano</i>	<b>Filo guida</b>
LV	<i>Latvian/Latviešu</i>	<b>Vadītājstīga</b>
LT	<i>Lithuanian/Lietuvių k.</i>	<b>Kreipiamoji viela</b>
NO	<i>Norwegian/Norsk</i>	<b>Ledesonde</b>
PL	<i>Polish/Polski</i>	<b>Prowadnik</b>
PT	<i>Portuguese/Potuguês</i>	<b>Fio-guia</b>
RO	<i>Romanian/Română</i>	<b>Fir de ghidaj</b>
RU	<i>Russian/Русский</i>	<b>Проволочный проводник катетера</b>
SK	<i>Slovak/Slovenčina</i>	<b>Vodiaci drôt</b>
SL	<i>Slovenian/Slovenščina</i>	<b>Vodilna žica</b>
ES	<i>Spanish/Español</i>	<b>Alambre guía</b>
SV	<i>Swedish/Svensk</i>	<b>ledare</b>
TR	<i>Turkish/Türkçe</i>	<b>Kılavuz Tel</b>

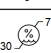
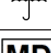
**EN: Explanation of symbols used on the package labels - BG: Значение на символите върху етикетите на опаковката - ZH: 包装标签上的标志说明 - HR: Objašnjenja simbola koji se koriste na naljepnici na pakiranju - CS: Vysvětlení symbolů použitých na štítcích obalů - DA: Forklaring af symboler anvendt på emballagens etiketter**

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## Accoat Guide Wire

### Description:

The device is an EO sterilized and non-pyrogenic guide wire.

Product name	Subtype	Tip subtype	Available variants
Accoat Guide Wire	Seldinger	<ul style="list-style-type: none"> <li>Straight</li> <li>J-tip</li> <li>Double ended</li> </ul>	<ul style="list-style-type: none"> <li>Stainless steel with 2 coating variants: with/without PTFE</li> </ul>
	Seldinger Movable	<ul style="list-style-type: none"> <li>Straight</li> <li>J-tip</li> </ul>	<ul style="list-style-type: none"> <li>Guide wire diameter</li> <li>Tip shape</li> <li>Guide wire stiffness</li> <li>Guide wire length</li> </ul>

For further information, see SP Medical's website.

### Intended purpose:

The Accoat Guide Wire is indicated for general intravascular and coronary arterial use to aid in the selective placement of interventional devices during diagnostic and / or therapeutic procedures.

### Contraindications:

#### Not for use in the cerebral vasculature.

It is always the physician's responsibility to determine and ensure the patient's suitability for the procedure where the Accoat Guide Wire is used.

### Warnings:

1. SINGLE USE. This device is intended for single use only. Reusing the device involves high risk of contamination and locking of the guide wire inside the interventional device or human body due to guide wire kinking.
2. This device is provided in sterile condition. Do not re-sterilize or re-use the device.
3. Do not withdraw PTFE coated guide wire through a metal cannula.
4. Never advance or withdraw the guide wire against resistance until the cause of the resistance is determined by fluoroscopy.
5. Do not attempt to move the guide wire without observing the resultant tip response.
6. To be used before the expiry date stated on the package.

### Precautions:

1. Prior to opening, the sterile package should be carefully inspected to see if it is still intact. The inspection must be carried out just before the medical device is used and in good light conditions. Sealing quality on the transparent side of the peel pouch must be carefully inspected for channels and cracks. Both front and back side are to be inspected for pinholes or other packaging defects such as abrasions or cracks. Use a magnification tool if in doubt of sterile barrier integrity.
2. Do not use if the package / sterile barrier is broken / damaged.
3. Prior to use, carefully inspect the guide wire for bends, kinks or other damages.
4. Do not use damaged product.
5. The device should only be used by experienced physicians, trained in invasive techniques, the use of guide wires and familiar with side effects and hazards commonly associated with interventional procedures.
6. When using movable-core guide wires, do not advance the movable core when the tip is in a curved shape.
7. Never twist or apply excessive force as the core might penetrate the coil and cause vessel damage.
8. To avoid arterial damage, make sure that the movable core is sufficiently retracted when the device is in the proximity of the heart.
9. The Accoat Guide Wire contains a metallic core, do not use with any inappropriate equipment (e.g., MRI).

### Potential side effects:

Possible complications include, but are not limited to the following:

- Vessel / tract trauma
- Vessel / tract narrowing / occlusion
- Infection and inflammation
- Allergic reaction

### Serious incident:

Every serious incident that has occurred in relation to the device should be reported without delay to the manufacturer and the competent authority of the Member State in which the user and / or patient is established. For further information please contact the manufacturer.

### Compatibility:

The Accoat Guide Wire may be used in combination with interventional devices in CathLab environment. The end-user is solely responsible for selection of appropriate device. There should be at least 0.0004" (0,01 mm) clearance between the lumen of the interventional device and the guide wire.

### Preparations for use:

1. Carefully remove the guide wire from the dispenser.
2. Inspect the guide wire thoroughly to make certain that it is not kinked or otherwise damaged.
3. Confirm the compatibility of the guide wire diameter with the interventional device before actual use.

### Directions for use:

1. When introducing the guide wire into the interventional device ensure that at least 2 centimeters of guide wire extend from the proximal hub. This will prevent the guide wire from slipping inside e.g., the catheter.
2. To aid in the selective placement of the interventional device into a particular vessel, gently rotate the proximal end of the guide wire as it is advanced forward.
3. To prevent contrast agent crystallization / clotting, a continuous saline flush should be maintained between the interventional device and the guide wire during the procedure. The size of the syringe used to flush the lumen should be adapted to the length and diameter of the interventional device.
4. Between uses, during the same procedure, place the guide wire in a saline-filled container, or fill the dispenser with saline, using the luer lock provided in the package and reinsert the guide wire in the dispenser, distal end first. Make sure to leave a segment of the proximal end outside the dispenser.

### Storage:

To be stored under cool, dark and dry conditions. To be kept away from heat. Humidity limits: 30% - 70%. Temperature limits: 2°C - 49°C

### Expiry date:

The Accoat Guide Wire must be used no later than stated on the label.

### Disposal:

Disposal shall be handled in accordance with accepted medical practice and applicable local laws and regulations for hazardous waste.

### Safety information:

The IFU is delivered in hard copy together with products and kept updated and available for electronic download on SP Medical's website. Summary of safety and clinical performance can be obtained via Eudamed using SP Medical as legal manufacturer and the basic UDI-DI: 5710086-GWACJU. <https://ec.europa.eu/tools/eudamed/#/screen/search-device>

### Disclaimer of liability:

SP Medical A/S is not liable for defects / deterioration resulting from abnormal use or modifications made to the product and under these circumstances not covered by the guarantee. SP Medical A/S disclaims liability for direct or indirect injuries that may occur as a consequence of the product being modified or wrongly used.

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